

Medical Disposable Products



Product Brochure

ABOUT US



JETBIO International Pvt. Ltd. is an Indian medical device manufacturing and export company dedicated to supplying high-quality, safe, and cost-effective healthcare products to global markets. With a strong commitment to quality, innovation, and customer satisfaction, we aim to support healthcare providers with reliable medical and surgical disposable products.

We specialize in the manufacturing and supply of hospital consumables, including IV cannula, Infusion sets, Haeomodialysis Blood Tubing Line, Yankauer suction set, Umbilical cord clamp, and other disposable medical products. Our production facilities follow strict quality control systems and international manufacturing standards to ensure consistent product performance and patient safety.

At JETBIO International, we focus on continuous improvement in product design, material selection, and packaging to meet the evolving needs of international healthcare markets. Along with competitive pricing, we provide efficient production planning, timely deliveries, and complete export documentation support.

We believe in building long-term partnerships with distributors, hospitals, and healthcare institutions worldwide by offering dependable products, responsive communication, and flexible OEM / private labeling solutions.

Our Mission

To deliver high-quality medical disposable products at competitive prices, making healthcare more accessible worldwide.

Our Vision

To become a globally trusted brand in medical and surgical disposable products.

INDEX

HYPODERMIC DISPOSABLE SYRINGE

IV. CANNULA / CATHETER

3 WAY STOP COCK

3 WAY STOP COCK WITH EXTENSION TUBE

EXTENSION TUBE

PRESSURE MONITORING LINE

INFUSION SET

BLOOD TRANSFUSION SET

MEASUREMENT VOLUME SET

IV FLOW REGULATOR

HEMODIALYSIS BLOOD TUBING LINE

AV FISTULA

SCALP VEIN SET

YANKAUER SUCTION SET

SUCTION CATHETER

UMBILICAL CORD CLAMP



JB VAN Hypodermic Disposable Syringe

Product Code: JB1401

- Single use only, E.T.O. Sterile, Non-Toxic, Non-Pyrogenic.
- Available in 1ml, 2ml, 3ml, 5ml, 10ml, 20ml & 50ml with & without needle.
- Available in Luer lock.
- Manufactured with Non-toxic medical grade Polypropylene and High Density Poly ethylene.
- Packed : Ribbon or Blister Package

Technical Specifications

Syringe	Needle Size	Packing Type	Packing Inner
1ml	26G×½", 26G×1½"	Ribbon / Blister	100 Pcs.
2ml	W/o Needle	Ribbon / Blister	100 Pcs.
2ml	21G, 22G, 23G, 24G, 26G	Ribbon / Blister	100 Pcs.
3ml	W/o Needle	Ribbon / Blister	100 Pcs.
3ml	21G, 22G, 23G, 24G, 26G	Ribbon / Blister	100 Pcs.
5ml	W/o Needle	Ribbon / Blister	100 Pcs.
5ml	21G, 22G, 23G, 24G, 26G	Ribbon / Blister	100 Pcs.
10ml	W/o Needle	Ribbon / Blister	100 Pcs.
10ml	21G, 22G, 23G	Ribbon / Blister	100 Pcs.
20ml	W/o Needle	Ribbon / Blister	50 Pcs.
20ml	20G, 21G, 22G	Ribbon / Blister	50 Pcs.
50ml	W/o Needle	Ribbon	25 Pcs.
50ml	20G, 21G, 22G	Ribbon	25 Pcs.



14G 16G 18G 20G 22G 24G 26G

JB FLON

I.V. Cannula / Catheter
with Wings & Injection Port

Product Code: JB0101

- Stainless steel, silicone tipped needle for smooth penetration.
- Injection port with one way silicone injection valve.
- Specially designed protection cap with a recessed plug to cover the injection port and minimize contamination.
- Angled and grooved wings for secure fixation.
- Radio-opaque lines for accurate radiographic detection.
- Luer cap for blocking the catheter when not in use.
- 6% luer taper for compatibility with all standard devices.
- Disposable, sterile and non-pyrogenic.



24G 26G

JB NEO

I.V. Cannula / Catheter
with Small wings & without Injection Port

Product Code: JB0102

- Small wings and absence of injection port permits easy handling.
- Large internal diameter for increased flow rate and reduced peel back effect.
- Siliconized needle tip for atraumatic tissue penetration.
- Transparent flashback chamber for easy visualization of blood.
- Available in size 24G & 26G.
- Disposable, sterile and non-pyrogenic.
- Sealed blister packing.

Technical Specifications

Gauge	Colour	IN. Ext. (in mm)	Length (in mm)	Water flow (in mm)
14	Orange	1.6-2.1	45	270
16	Grey	1.3-1.7	45	180
17	White	1.2-1.5	45	125
18	Green	0.9-1.3	45	90
18	Green	0.9-1.3	33	90
20	Pink	0.8-1.1	32	60
22	Blue	0.6-0.9	25	36
24	Yellow	0.5-0.7	19	20
26	Violet	0.4-0.6	19	13



JB Flex^{pro}

3 Way Stop cock with extension line

Product Code: JB0203

- Minimizes chances of mechanical irritation and infection by taking.
- The administration site away from insertion site.
- The integrated three-way stopcock offers multiple infusion lines.
- Smooth internal surface to minimize turbulence.
- Disposable, sterile and non-pyrogenic.

Technical Specifications

Length (cm.)	Int.- Ext. Dia. (mm.)	Qty. in carton Inner/ Outer
5	3.3-4.5	30/180
7	3.3-4.5	30/180
10	3.3-4.5	30/180
25	3.3-4.5	30/180
50	3.3-4.5	30/180
75	3.3-4.5	30/180
100	3.3-4.5	30/180
125	3.3-4.5	30/180
150	3.3-4.5	30/180
200	3.3-4.5	30/180
250	3.3-4.5	30/180

JB Flex

3 Way Stop cock

Product Code: JB0201 - JB0202

- Multiple channels for multiple infusion therapy.
- Designed to with stand pressure upto 5 bars.
- Smooth, fully rotatable tap.
- Stability tested pneumatically and hydrostatically.
- Tap turns every 360° without limitation.
- Flow shuts off every 90°.
- Minimal dead space in ports to ensure precise drug administration and maximum infusion flow. Continuous flow channels.
- Arrow indication marks on top to indicate direction of flow.
- Low profile for stable and safe positioning.
- Available in lipid resistant/non-lipid resistant varieties.
- Available with blue and red pegs for easier identification.
- Disposable, sterile and non-pyrogenic. Sealed blister packing.

Technical Specifications

Colour Code	Ref. No.		Qty. in carton Inner/ Outer
	Non- Lipid Resistant	Lipid Resistant	
Blue	JB0201BL	JB0202BL	50/500
White	JB0201WH	JB0202WH	50/500
Red	JB0201RE	JB0202RE	50/500



JB Ex

Extension Line

Product Code: JB0902

JB Ex^{pro}

Pressure Monitoring Line

Product Code: JB0901

- Made from medical grade PVC.
- Minimizes the movement at the cannulation site by taking the fluid administration site away from the insertion site.
- Smooth internal surface to minimize turbulence.
- Disposable, sterile and non-pyrogenic.

- Manufactured to withstand pressures upto 40 kg/cm².
- Fitted with male and female luer locks for secure connection.
- Smooth internal surface to minimize turbulence.
- Disposable, sterile and non-pyrogenic.

Technical Specifications

Length (cm.)	Int.- Ext. Dia. (mm.)	Qty. in carton Inner/ Outer
5	2.8-4.0	60/360
7	2.8-4.0	60/360
10	2.8-4.0	60/360
25	2.8-4.0	60/360
50	2.8-4.0	60/360
75	2.8-4.0	60/360
100	2.8-4.0	60/360
125	2.8-4.0	60/360
150	2.8-4.0	60/360
200	2.8-4.0	50/300
250	2.8-4.0	50/300

Technical Specifications

Int.- Ext. Dia. (mm.)	Length (cm.)	Qty. in carton Inner/ Outer
1.0-3.0	5	60/360
1.0-3.0	7	60/360
1.0-3.0	10	60/360
1.0-3.0	25	60/360
1.0-3.0	50	60/360
1.0-3.0	75	60/360
1.0-3.0	100	60/360
1.0-3.0	125	60/360
1.0-3.0	150	60/360
1.0-3.0	200	50/300
1.0-3.0	250	50/300




JB Line

Hemodialysis Blood Tubing Line

Product Code: JB0801-05

- Blood tubing set has been specially designed for the patient to connect external system which extracts patient's blood to the dialyzer and then returns patient's blood from the dialyzer.
- Consists of two parts - arterial and venous - used during dialysis and attached with fistula and dialyzer.
- Unique chambers reduce foaming, increase air removal and do not trap epo compliant material of transducer protectors. This helps in ensuring secure machine fit and less incidence of wet out.
- Color coded ergonomic clamps and dialyzer connectors with tethered caps that seal tightly for disposal.
- Kink resistant blood tubing handles easily and ensures fast priming
- Additional standard accessories
- Sterilization : ETO
- Individual Pouch Pack, Outer: 50 Pcs Pack

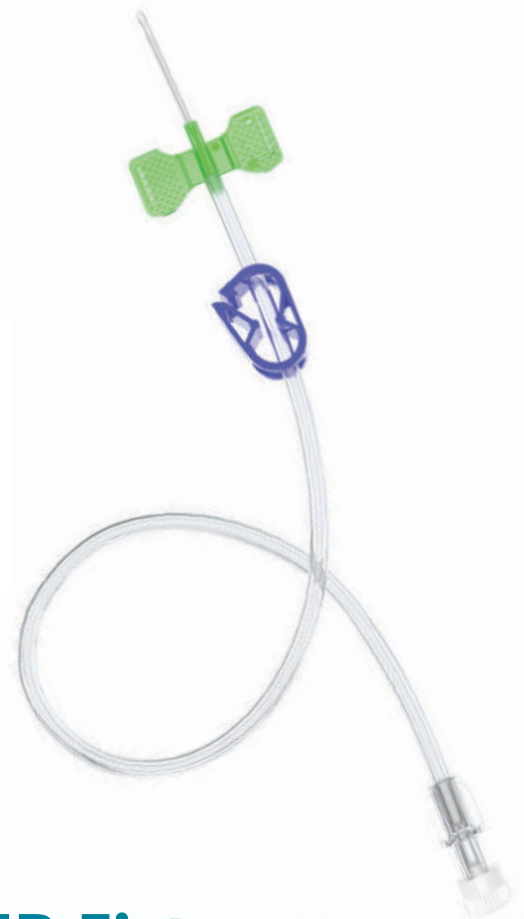

JB Fist

AV Fistula Needle

Product Code: JB1101-03

- Ultra-siliconized thin-walled sharp needle with bevelled back eye for easy insertion and minimize trauma to the patient
- Innovative design ensures patient safety and user comfort.
- Soft and kink-resistant medical-grade PVC tubing.
- Color-coded flexible wings for easy identification. Equipped with a 6% female luer taper port for a secure connection with other medical devices, preventing any leakage.
- Needle Gauge 15G, 16G and 17G
- With or Without Back-eye Needle
- Twin Set (Arterial & Venous)
- Sterilization : Gamma Irradiation or ETO
- Pack Size : Inner : 50 Sets Pack,
Outer : 500 Sets Pack



JB Set

Infusion Set - V / NV
with Latex bulb / with latex tube
Product Code: JB0301



JB Set^{pro}

Infusion Set - V
with Luer Lock & Y-Injection port
Product Code: JB0302



JB Set⁺

Infusion Set - V
with Luer Lock & Y-Injection port
with auto stop chamber
Product Code: JB0303



JB Set^{Blood}

Blood Transfusion Set
Product Code: JB0304



JB Set^{BURETTE}
Measured Volume Set
Product Code: JB0305

For Infusion Set

- Transparent & flexible drip chamber
- 15 micron fluid filter in drip chamber(if required)
- Soft and kink resistant PVC tubing
- Smooth roller clamp facilitates easy & controlled adjustment of fluid rates
- Standard tube length : 150 cms

For BT Set

- For transfusion of blood or blood components
- Transparent & flexible drip chamber
- Approximately 20 drops/ml
- Soft and kink resistant PVC tubing
- With 200 micron blood filter
- Also available in Double drip chamber
- With or without needle

JETBIO

Yankauer Suction Set
Product Code: JB0401



JETBIO

Suction Catheter
Product Code: JB0501



JETBIO

Scalp Vein Set
Product Code: JB1001



JB FIO

IV Flow Regulator
Product Code: JB1301



JETBIO

Umbilical Cord Clamp
Product Code: JB0601

NOTES

This image shows a full page of white paper with horizontal blue ruling lines. The lines are evenly spaced and run across the width of the page. There is no handwriting or other markings on the paper.

Your Health Our Priority



JETBIO International Private Limited

Plot No. 172, Second Floor PKT K Sec-3, Bawana DISDC Landmark

Bawana Industrial Complex Delhi 110039

Customer Care No. + 91 9625550899

Email : info@jetbio.in, Web.: www.jetbio.in

